

Chairside Attachment Pickup Material



Chairside Attachment Pickup Material Kits. Each kit contains:



Syringe and Plunger (8g)

15 Mixing Tips

15 Angled Tips

Features

Optimal handling characteristics

•Tissue-colored material bonds to itself

•Packaged with a straight and angled tip for precise application

•Dual cure - self or light cure

•The material will set in approximately 5 minutes from the moment the material is expressed from the tip. The working time is approximately 2 minutes and 15 seconds.

•The material will self-cure in approximately 7 minutes from the time the material is expressed from the tip. We recommend waiting the full 7 minutes to ensure reliable attachment fixation. Curing time can be shortened by light curing.

Cost efficient


•Shortened mixing tip length reduces waste and cost and increases the number of Attachments processed per syringe

•No primer required, reducing associated chairside time and cost

Designed with patient comfort in mind

•Odorless and tasteless •Low cure temperature

Six Easy Steps

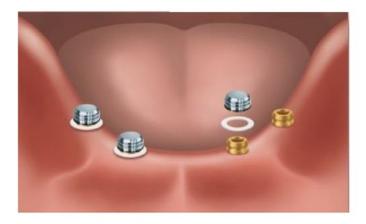

1. PLACE

Place blockout spacers around and firmly seat retention caps on to abutments.




2. APPLY

Apply fit check marking/indicating paste inside of the denture. Insert the denture (in the mouth) in position over the Retention Caps. Mark the areas where the denture will need to be relieved to allow space for the caps to be picked up



3. RELIEVE

Relieve the areas marked with an acrylic bur. Try in the denture to verify that the Retention Caps are not in contact with the acrylic in any area. <u>Cut lingual/palatal</u> <u>escape vent windows in the denture</u> to visualize full seating and for an excess material vent. Cut an <u>undercut</u> around the circumference of the relief areas for mechanical retention

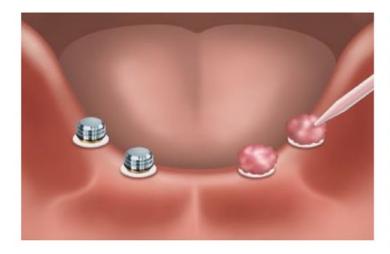

4. DRY

Dry the Retention Caps. Apply a small amount of CHAIRSIDE Material around the circumference of eacl cap. Place CHAIRSIDE Material into the relief areas of the denture and seat the denture over the caps and onto the tissue. Have the patient close into <u>light occlusion</u>--or hold the denture in place with finger pressure over the attachments-- and hold while the CHAIRSIDE Material sets. The set time of the material from the time the material is mixed to the time it is set i five (5) minutes. The working time for the material is two (2) minutes, fifteen (15) seconds. Light cure to expedite setting.



5. DISENGAGE

Disengage the denture from the attachments and removies from the mouth. Verify that the Retention Caps have been securely processed into the denture. Fill any voids light cure and polish the denture. The material will bond to itself and will cure within 30 seconds with light application.



6. FIT

Reseat the denture in the mouth and press down to engage the males on the Attachments.